




technoclone

THROMBIN GENERATION

Fully automated - as easy as a routine assay

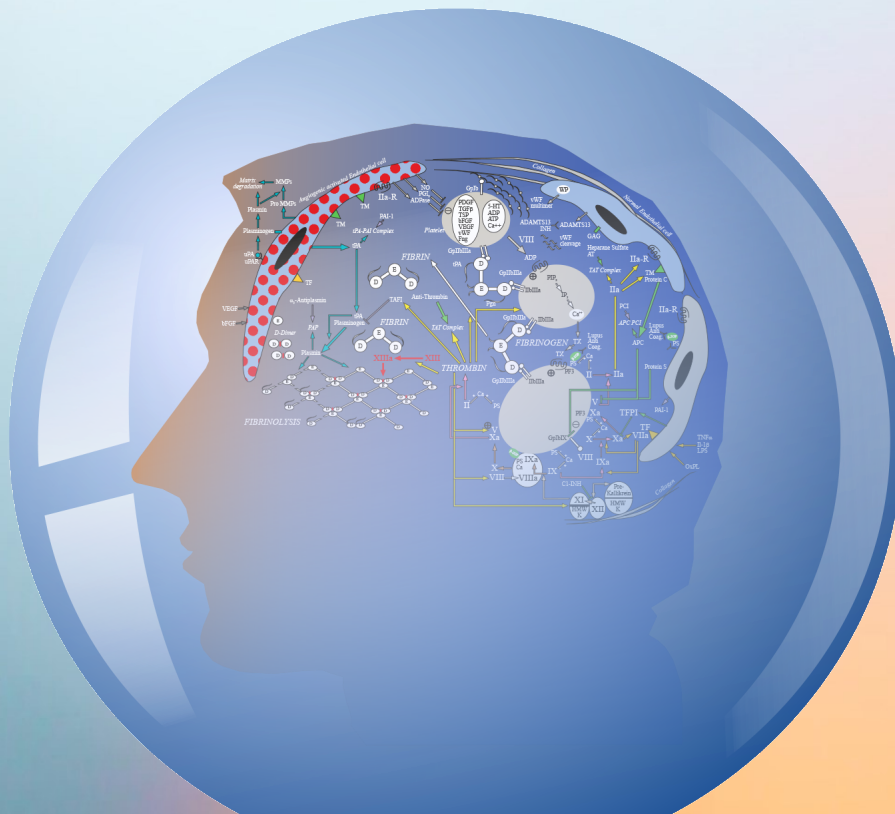

CE marked

A clear view into the future!

Thrombin Generation on Ceveron® t100 / Ceveron® s100

The **Ceveron® TGA Kits**, for the measurement of **bleeding tendencies**, **thrombophilic tendencies** or **anticoagulated patients**, are optimized for the fully automated **CE marked** measurement of Thrombin Generation on **Ceveron® t100** and **Ceveron® s100**:

- **Easy** as a routine assay
(Measurement from the same patient sample as a routine parameter)
- **Parallel testing** of routine & thrombin generation samples
- **Precise** (intra/inter-assay CV < 5%)
- **Fast** (short assay time of ~20 min for Peak Thrombin)
- Automated **normalisation**
- Up to **36 samples** in a single run

Product	Description	REF	Package
Ceveron® TGA RB Kit	For measurement of bleeding tendencies. Low concentration of phospholipid micelles containing low rhTF. the kit contains: 2 x 1 ml Ceveron® TGA Reagent RB 2 x 3 ml Ceveron® TGA Substrate 2 x 1.5 ml Ceveron® TGA Reaction Buffer 2 x 1 ml Ceveron® TGA Control Low 2 x 1 ml Ceveron® TGA Control High 1 x 25 ml CaCl₂ 25 mM	5006011	120 tests
Ceveron® TGA RC Low Kit	For measurement of thrombophilic tendencies. Low concentration of phospholipid micelles containing rhTF. the kit contains: 2 x 1 ml Ceveron® TGA Reagent RC Low 2 x 3 ml Ceveron® TGA Substrate 2 x 1.5 ml Ceveron® TGA Reaction Buffer 2 x 1 ml Ceveron® TGA Control Low 2 x 1 ml Ceveron® TGA Control High 1 x 25 ml CaCl₂ 25 mM	5006013	120 tests
Ceveron® TGA RC High Kit	For measurement in anticoagulated patients. High concentration of phospholipid micelles containing rhTF. the kit contains: 2 x 1 ml Ceveron® TGA Reagent RC High 2 x 3 ml Ceveron® TGA Substrate 2 x 1.5 ml Ceveron® TGA Reaction Buffer 2 x 1 ml Ceveron® TGA Control Low 2 x 1 ml Ceveron® TGA Control High 1 x 25 ml CaCl₂ 25 mM	5006015	120 tests
Calibrator Set			
Ceveron® TGA CAL Set	1 x 3 ml Ceveron® TGA Buffer	5006347	1 Set

